

# Gomco® 405 Aspirator



Model 01-12-4005 comes with one 1500 ml disposable canister, suction tubing, and 3 filters

The Gomco AC Model 405 is a high-performance general use tabletop aspirator for use in hospitals, clinics, physician's offices, nursing home facilities and home care. This durable unit meets the need for controlled vacuum regulation up to 635 mm Hg of vacuum with a flow rate of 40 lpm at open flow. The aspirator has a high-impact plastic housing and powerful rotary vane pump, reduced operating noise level, and weighs less than 15 pounds. Each unit comes with either a 1500 ml disposable canister or a 1200 ml reusable glass bottle, a suction tubing kit, and three hydrophobic bacteria filters.


## Standard Features

- Powerful rotary vane pump excels at pharyngeal and nasopharyngeal procedures
- Compact ... the 405 takes up less than one square foot of space
- Self-contained unit provides protection against pump and gauge contamination, and allows for easy clean-up
- Low maintenance requirements
- Durable high-impact plastic case
- Disposable hydrophobic bacteria filter eliminates pump flooding and helps prevent fluid and aerosol contamination

# Gomco® 405 Aspirator



01-90-3100

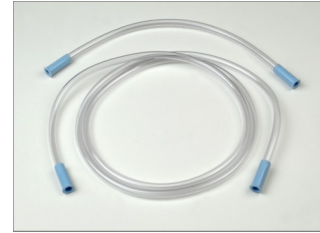

S610100



20-08-0003



01-90-3489

Model 01-12-0405 comes with one 1200 ml glass collection bottle, suction tubing, and 3 filters

## SPECIFICATIONS

	<u>01-12 Models</u>	<u>01-23 Model</u>
<b>Vacuum Range:</b>	-80 to -635 mm Hg	-25 to -555 mm Hg
<b>Flow Rate:</b>	40 lpm at open flow	25 lpm at open flow
<b>Sound Level:</b>	Less than 54 dBA	
<b>Dimensions:</b>	12.5" W x 10" D x 9.5" H (31.75 x 25.4 x 24.13 cm)	
<b>Weight:</b>	14.5 lbs. (7 kg)	
<b>UL/ETL Listed:</b>	Yes	
<b>Warranty:</b>	1 year	
<b>Regulatory Approvals:</b>	FDA, CE, HC	

## ORDERING INFORMATION

<b>01-12-0405</b>	Aspirator, 1200 ml glass bottle, suction tubing kit, 3 filters, 115V/60 Hz
<b>01-12-4005</b>	Aspirator, 1500 ml disposable canister, suction tubing kit, 3 filters, 115V/60 Hz
<b>01-23-0405</b>	Aspirator, 1200 ml glass bottle, suction tubing kit, 3 filters, 230V/50-60 Hz

## ACCESSORIES

<b>20-08-0003</b>	Disposable Collection Canister, 1500 ml, stem inlet (48/case)
<b>20-08-0016</b>	Disposable Collection Canister, 1500 ml, stem inlet (16/case)
<b>01-90-3489</b>	Glass Collection Bottle, with metal cap & float assembly, 1200 ml
<b>01-90-3497</b>	Glass Collection Bottle (bottle only), 1200 ml (6/case)
<b>01-90-3487</b>	Metal Cap & Float Assembly for 1200 ml glass bottle
<b>01-90-3100</b>	Disposable Hydrophobic Bacteria Filter, 1/8" NPT x 1/4" hose barb (3/pkg.)
<b>S610100</b>	Disposable Suction Tubing Kit: 1 each 15" & 72", 1/4" ID, blue-tipped, clear PVC



Allied Healthcare Products, Inc. ▪ 1720 Sublette Avenue ▪ St. Louis, MO 63110 USA  
314/771-2400 ▪ 800/444-3954 ▪ FAX: 314/771-0650 ▪ www.alliedhpi.com

U.S.  
Toll Free: 800/444-3940  
FAX: 800/477-7701

Canada & Mexico  
Toll Free: 800/446-0552  
FAX: 800/246-6201

International  
Telephone: 314/268-1683  
FAX: 314/771-5183

ISO 13485: 2003

©Allied Healthcare Products, Inc.®

All specifications are nominal and subject to change without notice.

Form 21-00-0004  
January 2014