


# Constant Flow Selector Valves

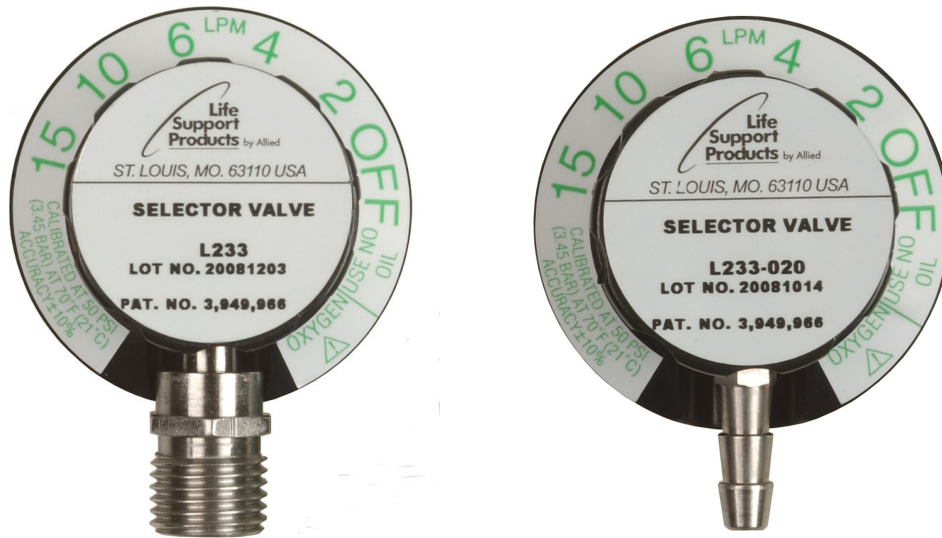

LSP Constant Flow Selector Valves are designed for use with hospital and/or portable oxygen systems where a metered supply of oxygen is required.


Unlike standard flow meters and flow gauges, the LSP selector valves have five distinct flow rate settings. These highly accurate selector valves offer repeatable flow rates regardless of orientation or visibility.

The unit has 1/8" NPT thread to allow attachments of Schrader, Ohio Medical, Hansen or NCG disconnects. Constant flow selector valves lock into place with each selected liter flow setting. All knobs are clearly marked and are accurately adjustable, even in low light situations.

The polycarbonate selector knob is virtually unbreakable and designed to withstand rugged abuse. Lightweight aluminum and a low-profile design make these valves ideal for emergency medical applications.

The selector valve's ease of operation in demanding environments and rugged low-profile design have made it the choice of emergency medical professionals.

# Constant Flow Selector Valves

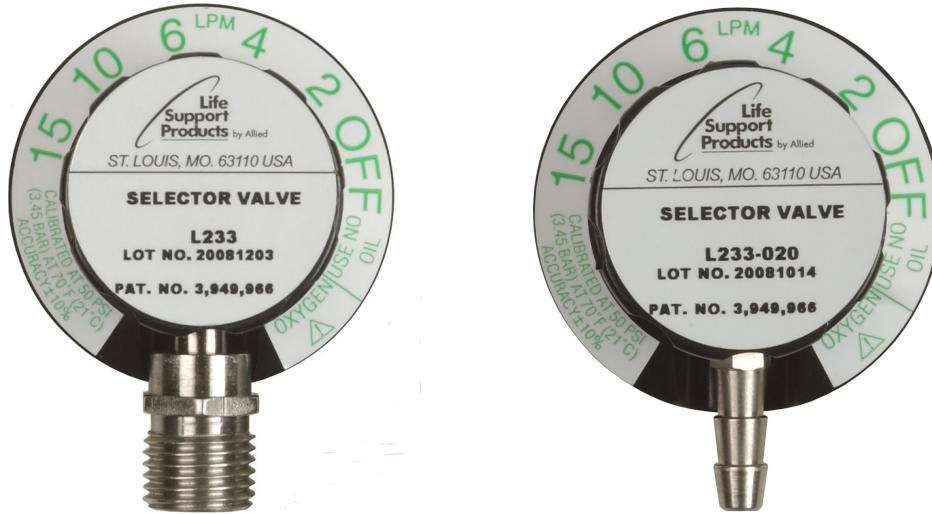

## SPECIFICATIONS

- Inlet Pressure:** 50 psig nominal  
**Operating Temperature:** -30° F to 125° F  
**Storage Temperature:** -40° F to 160° F  
**Flow Rates:** 2, 4, 6, 10 and 15 lpm (all flow rate settings accurate within  $\pm 10$  per cent)  
**Filter:** 25 micron stainless steel mesh screen  
**Inlet Port:** 1/8" NPT female  
**Outlet Fittings:**  
    **L233 & L233-050:** Barbed-type fitting which can be removed for access to a standard male oxygen DISS-type fitting  
    **L233-020:** Barbed-type fitting only  
**Material:** Body - anodized aluminum  
    Knob - polycarbonate

## ORDERING INFORMATION

- L233** Constant Flow Selector Valve with oxygen DISS outlet, 2, 4, 6, 10 & 15 lpm  
**L233-020** Constant Flow Selector Valve with oxygen barbed outlet, 2, 4, 6, 10 & 15 lpm  
**L233-050** Constant Flow Selector Valve, with oxygen DISS outlet, 2, 4, 6, 10 & 25 lpm



**EMERGO EUROPE**  
Molenstraat 15  
2513 BH, The Hague  
The Netherlands



Allied Healthcare Products, Inc. • 1720 Sublette Avenue • St. Louis, MO 63110 USA  
314/771-2400 • 800/444-3954 • FAX: 314/771-0650 • www.alliedhpi.com

U.S.  
Toll Free: 800/444-3940  
FAX: 800/477-7701

Canada & Mexico  
Toll Free: 800/446-0552  
FAX: 800/246-6201

International  
Telephone: 314/268-1683  
FAX: 314/771-5183

ISO 13485: 2003

©Allied Healthcare Products, Inc.®

All specifications are nominal and subject to change without notice.

Form 90-00-0018  
Rev. January 2014