

Getting Good Connections Has Never Been Easier.




Connect2® Medical Gas Outlet Station: The Right Connection for New or Existing Structures



- The only medical gas outlet that accepts both Chemetron and Ohmeda adapters (patent #5,562,121).
- Push-button mechanism for quick release of Chemetron adapter; Ohmeda adapter has a twist mechanism for quick release.
- Chemetron adapter can support secondary equipment, including a suction regulator with a half-gallon collection bottle.
- Eliminates the need for hospital personnel to carry both styles of adapters to ensure proper connection.
- Fail-safe color coding and safety keyed to ensure appropriate gas.
- Valve body O-ring easily accessible for maintenance; inlet tube has secondary check unit to permit removal of valve body without shutting off gas to the outlet.
- Cleaned for oxygen use.
- UL and CSA* listed. NFPA compliant.

Ordering Information		
Connect2®	Domestic Quick Connect	ISO* Quick Connect
Oxygen	64-02-5001P	64-02-5051P
Vacuum	64-02-5002P	64-02-5052P
Air	64-02-5003P	64-02-5053P
Nitrous Oxide	64-02-5004P	64-02-5004P
WAGD	64-02-5006P	64-02-5006P
Carbon Dioxide	64-02-5007P	64-02-5007P
O2/CO2<7%	64-02-5008P	64-02-5058P

*CSA listed with ISO colors only.

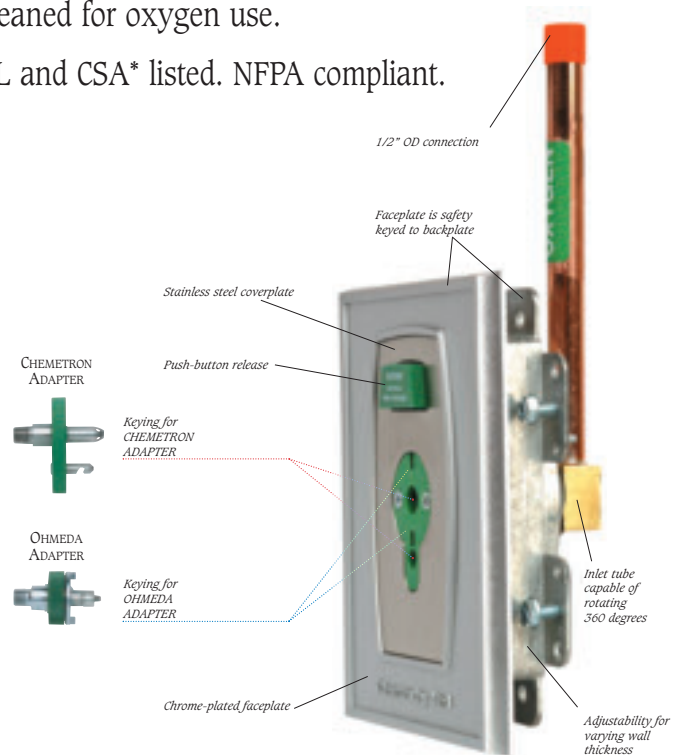

The Chemetron® Outlet: An Industry Standard for More Than 40 Years



- Release mechanism designed for one-handed thrust mounting and fingertip release of secondary equipment.
- Prohibits the interchange of key components, backbox to adapters.
- Color coded for error-proof identification.
- Adjustable modular design backbox with rotating inlet tube for varying wall thickness and easy installation.
- Conforms with NFPA and CGA standards.
- Cleaned for oxygen use.
- UL and CSA listed. NFPA compliant.
- Fingertip push-button release mechanism for

Ordering Information				
Chemetron®	Domestic		ISO	
	Quick Connect	DISS	Quick Connect	DISS
Oxygen	64-01-5001P	64-04-5601P	64-01-5051P	64-04-5651P
Vacuum	64-01-5002P	64-04-5602P	64-01-5052P	64-04-5652P
Air	64-01-5003P	64-04-5603P	64-01-5053P	64-04-5653P
Nitrous Oxide	64-01-5004P	64-04-5604P	64-01-5004P	64-04-5604P
WAGD	64-01-5006P	64-04-5606P	64-01-5006P	64-04-5606P
Carbon Dioxide	64-01-5007P	64-04-5607P	64-01-5007P	64-04-5607P
O2/CO2<7%	64-01-5008P	64-04-5608P	64-01-5058P	64-04-5658P
Nitrogen	N/A		N/A	

The Oxequip® Outlet: Flexibility That Meets All the Standards

- Outlet design meets domestic and international standards for medical gas equipment.
- Modular design gives installer flexibility to vary sequence of gases for individual rooms.
- Color-coded outlet has non-interchangeable indexing system to prevent accidental interchange.
- Improved secondary check increases gas flow to patient.
- Easier installation due to adjustability for varying wall thickness and adjustable rotating inlet tube.
- Chrome-plated die-cast faceplate provides durability and improved appearance.
- Cleaned for oxygen use.
- UL and CSA listed. NFPA compliant.

Ordering Information				
Oxequip®	Domestic		ISO	
	Quick Connect	DISS	Quick Connect	DISS
Oxygen	1930-2111	1930-1111	1930-2121	1930-1121
Vacuum	1930-2211	1930-1211	1930-2221	1930-122
Air	1930-2411	1930-1411	1930-2421	1930-1421
Nitrous Oxide	1930-2311	1930-1311	1930-2321	1930-1321
WAGD	1930-2811	1930-1911	1930-2821	1930-1921
Carbon Dioxide	1930-2511	1930-1511	1930-2521	1930-1521
Nitrogen	N/A	1930-1611	N/A	1930-1621



The PB Style Outlet: Easy to Use and Reliable



secondary equipment.

- Utilizes Chemetron backbox for easy upgrade to Chemetron completion kit.
- Die-cast faceplate is chrome-plated for long life and durability.
- Non-interchangeable indexing system prevents mislabeling.
- Meets NFPA and CGA standards.
- Cleaned for oxygen use.
- UL listed. NFPA compliant.

Ordering Information

	Domestic Quick Connect
PB Style	
Oxygen	64-03-5001
Vacuum	64-03-5002
Air	64-03-5003
Nitrous Oxide	64-03-5004



Allied Healthcare Products, Inc.™ • 1720 Sublette Avenue • St. Louis, MO 63110
314/771-2400 • 800/444-3954 • Fax 314/771-0650 • www.alliedhpi.com
International offices in Europe, Asia and South America

U.S.

Toll Free: 800/444-3960
Fax: 314/771-7537

Canada & Mexico

Toll Free: 800/446-0552
Fax: 800/246-6201

International

Telephone: 314/268-1683
Fax: 314/771-5183